


SOP for FS-Avoidance behavior chambers built by Zengyou

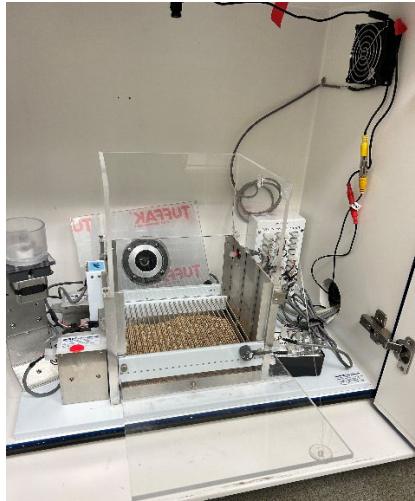

1. Environmental Habituation

- Bring the mice to the testing room and allow them to habituate to the environment for at least **30 minutes**.

2. Turn On Equipment

- Turn on the computer and cameras using the power switch.

3. Set Up Camera Monitoring

- Open **Blue Iris** and confirm the camera views for **Boxes 5–8**.

4. Set Up Med-PC

- Open **Med-PC**.
- Enter/check the initial parameters for each cage, including the **subject ID** and **protocol**.

5. Pre-Test Check

Before starting, confirm:

- **Syringe pump**: Ensure **0.2% saccharin solution** is loaded and all tubing is filled.
- **FS device**: Confirm the setting is approximately **0.3 mA**.

6. Start Behavioral Testing

- Place the mice in the behavioral chambers.
- Start the appropriate Med-PC program for each mouse.

7. After Testing

- Return the mice to their home cages.
- Record the **number of rewards** and **FS count** for each mouse.
- Measure and record the **remaining volume of saccharin solution**.

8. Clean the Chambers

- Clean the chambers with **70% ethanol**.
- Allow the chambers to dry before starting the next batch.